

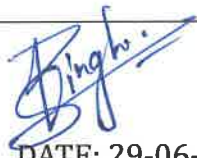
INSTRUCTIONS FOR USE

SPMM/RA/IND/43/ENFS/IFU, Rev. 00. Date: 29/06/2026



INSTRUCTIONS FOR USE

Enteral Feeding (Enfit) Syringe

 DATE: 29-06-2026	 DATE: 29-06-2026	 DATE: 29-06-2026
Prepared By Executive	Reviewed By Head	Approved By Director/Auth. Signatory

INSTRUCTIONS FOR USE

SPMM/RA/IND/43/ENFS/IFU, Rev. 00. Date: 29/06/2026

Device Description-

The Enteral Feeding (ENFit) Syringe is a single-use, three-piece syringe comprising a transparent graduated barrel, plunger rod, elastomeric plunger gasket, ENFit connector and protective tip cap. The threaded ENFit connector provides a secure connection to compatible enteral feeding tubes, medication ports and accessories.

The transparent barrel and clearly printed graduation scale support visual inspection and accurate measurement. The purple enteral identification colour helps distinguish the device from syringes intended for other clinical routes. Smooth plunger movement enables controlled aspiration and administration of fluids.

Sterilization: Non-Sterile / EO Sterile

Class of Medical Device: Class B



Fig- Enteral Feeding (Enfit) Syringe

Models/Variants-

Capacity	Sterility	Packing
1 mL	Non-Sterile / EO Sterile	Multi pack/individual Blister pack
3 mL	Non-Sterile / EO Sterile	Multi pack/individual Blister pack
5 mL	Non-Sterile / EO Sterile	Multi pack/individual Blister pack
10 mL	Non-Sterile / EO Sterile	Multi pack/individual Blister pack

INSTRUCTIONS FOR USE

SPMM/RA/IND/43/ENFS/IFU, Rev. 00. Date: 29/06/2026

20 mL	Non-Sterile / EO Sterile	Multi pack/individual Blister pack
30 mL	Non-Sterile / EO Sterile	Multi pack/individual Blister pack
50mL	Non-Sterile / EO Sterile	Multi pack/individual Blister pack
60 mL	Non-Sterile / EO Sterile	Multi pack/individual Blister pack
100mL	Non-Sterile / EO Sterile	Multi pack/individual Blister pack

Material-

- Medical-grade polypropylene (PP) - barrel and plunger rod
- Medical-grade polypropylene (PP) - ENFit connector and protective cap, as applicable
- Synthetic rubber / thermoplastic elastomer (TPE) - plunger gasket
- Medical-grade printing ink - graduation and identification markings
- Silicone lubricant - if applicable to the validated design

Intended Use

Enteral Feeding (Enfit) Syringe is used for the transport and administration of non-injectable liquid medication. The Tip Cap is used to cap the syringe to prevent leakage during transportation.

Intended User

- Pharmacists, Nurses and Physicians experienced in the use of the syringe.
- Experience levels will be from novice to expert.

Principle of Operation

Manual withdrawal of the plunger creates negative pressure within the barrel for aspiration of enteral fluid. Forward movement of the plunger generates positive pressure for controlled delivery through the ENFit connector. The threaded connector engages with a compatible ENFit male connector to provide a secure enteral connection.

INSTRUCTIONS FOR USE

SPMM/RA/IND/43/ENFS/IFU, Rev. 00. Date: 29/06/2026

Instructions For Use

1. Accurately measure prescribed dosage using the top edge of the stopper. Syringe scales are indicated in milliliters (mL).
2. After filling syringe, attach tip cap which is provided separately in small bag.
3. Before dispensing, REMOVE TIP CAP.
4. To avoid choking, slowly dispense the medication by pushing gently on the plunger rod.
5. Discard syringe properly after use.

Warnings

- Do not use if the product or package has visible particulates.
- Potential Procedure Related Complications: Exposure to toxic medication.

Precautions

- Pharmaceutical products prepared into these syringes should be administered as soon as possible. For light sensitive medication, use purple colored Enteral Feeding (Enfit) Syringe.
- Re-use may lead to infection or other illness/injury.

Shelf Life: 5 Years

Storage Conditions

Store in a clean, cool and dry place, protected from direct sunlight, excessive heat, moisture and contamination. Maintain the device in its original unopened packaging until use.

Disposal

Discard syringe properly after use according to applicable local waste disposal requirements.

EXPLANATION OF SYMBOLS-

INSTRUCTIONS FOR USE

SPMM/RA/IND/43/ENFS/IFU, Rev. 00. Date: 29/06/2026

Symbol	Explanation	Symbol	Explanation
	Caution		Catalogue Number
	Manufacturer		Batch code
	Date & Country of manufacturer		Use-by date
	Medical Device		Do not use hook
	Consult instructions for use		Do not use if package is damaged and consult <i>instructions for use</i>
	Keep away from sunlight		Recycling
	Keep dry		Handle with care
	This side up		Model number

Revision History

Rev. No.	Date	Description of Change
00	29/06/2026	Initial release of Instructions for Use for Enteral Feeding (Enfit) Syringe