



INSTRUCTIONS FOR USE

Nebulizer Mask Kit


DATE: 08/09/2026	DATE: 08/09/2026	DATE: 08/09/2026
Prepared By Executive	Reviewed By Head/MR	Approved By Director/Auth. Signatory

INSTRUCTIONS FOR USE

SPMM/RA/IND/70/NMK/IFU, Rev. No.: 00, Date: 08/09/2026

Product Name- Nebulizer Mask Kit

Image of the Device-

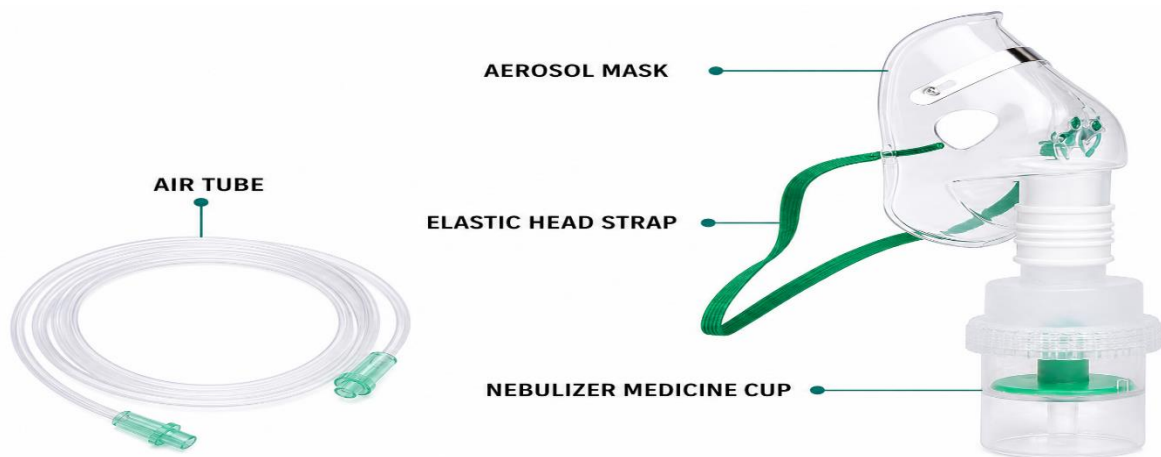

Intended Use-

The kit is intended to connect to a compatible nebulizer/compressor system to facilitate aerosol delivery of prescribed nebulized medication to a patient through an adult or pediatric face mask.

Warnings and Precautions

- Use only with compatible nebulizer equipment and according to the medication prescription.
- Do not use if the package is damaged, opened or contaminated. For sterile versions, do not use if the sterile barrier is compromised.
- Do not exceed the validated chamber fill volume.
- Ensure all connectors are securely fitted before starting therapy.
- Avoid kinking or occluding the connecting tube during use.
- Use the correct adult or pediatric mask size to obtain an appropriate facial fit.
- Dispose of used components in accordance with hospital policy and local regulations.
- This device is not waterproof or water resistant. Do not allow this device to become submerged in water.
- The device should not be touched with wet hands.
- Do not use the device in the presence of flammable anesthetic substances.

INSTRUCTIONS FOR USE

SPMM/RA/IND/70/NMK/IFU, Rev. No.: 00, Date: 08/09/2026

- Unplug the device before cleaning.
- Unplug the unit when it is not in use.
- Do not untwist the nebulizer cup during operation.
- Do not leave solution in the nebulizer cup and always empty when it is not in use.
- Place the device on a flat surface.
- This product is meant for intermittent use.
- In order to prevent the risk of cross contamination, the accessories should not be shared by different people.
- Children may only use the device under adult supervision.
- All accessories included are for Single Patient Use Only.

Storage Conditions-

Store in a clean, dry location, protected from direct sunlight, excessive heat, moisture and physical damage.

Directions for Use

- Open the package and inspect the kit for damage, contamination or loose connections.
- Add the prescribed volume of medication to the nebulizer chamber without exceeding the validated maximum fill volume.
- Close the chamber securely and connect the tube to the compatible nebulizer/compressor outlet.
- Attach the adult or pediatric mask to the chamber outlet.
- Place the mask comfortably over the patient's nose and mouth.
- Plug in and switch ON the nebulizer equipment to start treatment.
- Turn off the nebulizer equipment and unplug it once the treatment has finished.
- After therapy, dispose of used single-use components according to local requirements and applicable disposal procedures.

Shelf Life- 3 Years

EXPLANATION OF SYMBOLS:

Symbol	Explanation	Symbol	Explanation
	Caution		Catalogue Number

INSTRUCTIONS FOR USE

SPMM/RA/IND/70/NMK/IFU, Rev. No.: 00, Date: 08/09/2026

	Manufacturer		Batch code
	Date & Country of manufacturer		Use-by date
	Do not re-sterilize		Medical Device
	Consult instructions for use		Do not use if package is damaged
	Keep away from sunlight		Keep dry
	Do not re-use		Sterilized using Ethylene Oxide

Revision History:

Rev. No	Reason for Change	Date
00	Initial Release of IFU	08/09/2026