INSTRUCTION FOR USE

SURGICAL GOWN

INSTRUCTIONS FOR USE

SPMM/RA/IND/09/SRGG/IFU, Rev. 00, **Date:**26-09-2025

PRODUCT DESCRIPTION

The Surgical Gown is a sterile, single-use protective garment made from non-woven SMS (Spunbond–Meltblown–Spunbond) and polypropylene materials. It is designed to be worn by healthcare professionals during surgical procedures.

INTENDED USE / INTENDED PURPOSE

The Surgical Gown is intended for use during surgical procedures to help protect patients and healthcare personnel from the transfer of microorganisms, body fluids, and particulate materials.

INSTRUCTIONS FOR USE

1. **Inspect:** Check the packaging for damage and verify the expiry date before use.
2. **Open:** Perform hand hygiene and open the package using aseptic technique.
3. **Wear:** Unfold the gown and insert both arms into the sleeves without touching the outer surface.
4. **Secure:** Fasten the gown at the neck and waist with assistance, following standard sterile gowning procedures.
5. **Remove:** After use, carefully remove the gown by turning it inside out to avoid contact with the contaminated outer surface.
6. **Dispose:** Discard the used gown in accordance with applicable biomedical waste management regulations.

WARNINGS AND PRECAUTIONS

- For single use only. Do not reuse or resterilize.
- Do not use if the packaging is damaged or opened.
- Do not use after the expiry date.
- Do not use if the gown is torn, damaged, or visibly contaminated before use.
- Select an appropriate size to ensure proper coverage.
- Replace the gown if its protective integrity is compromised during use.

INSTRUCTIONS FOR USE

SPMM/RA/IND/09/SRGG/IFU, Rev. 00, **Date:**26-09-2025

STORAGE CONDITIONS

Store in a cool, dry, and clean place. Protect from direct sunlight, moisture, and physical damage.

SHELF LIFE

5 years from the date of manufacture, subject to the validated shelf life and expiry date indicated on the product label.

STERILIZATION

Sterilized using Ethylene Oxide (EO).

DISPOSAL

Dispose of the used gown in accordance with the local laws and regulations.

LEGAL MANUFACTURER AND MANUFACTURING FACILITY

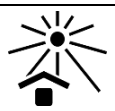
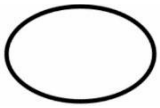
Manufacturer:
SPM Medicare Pvt, Ltd
B-40, Phase-II, Noida, Gautam Buddha Nagar, UP-201305, India.
Manufacturing Site 1- B-40, Phase-II, Noida,
Gautam Buddha Nagar, UP-201305, India.
Customer care: No.: +91-840-7070-718
Email: info@spmmedicare.com
Web Site: www.spmmedicare.com

EXPLANATION OF SYMBOLS:

Symbol	Explanation	Symbol	Explanation
	Caution		Catalogue Number
	Manufacturer		Batch code
	Date & Country of manufacturer		Use-by date
	Do not re-sterilize		Model number

INSTRUCTIONS FOR USE

SPMM/RA/IND/09/SRGG/IFU, Rev. 00, **Date:**26-09-2025

	Consult instructions for use		Do not use if package is damaged
	Keep away from sunlight		Non-Pyrogenic
	Do not re-use		Sterilized using Ethylene Oxide
	Keep dry		Handle with care
	Medical Device		Single Sterile Barrier System

Revision History:

Rev. No	Reason for Change	Date
00	Initial Release	26-09-2025